



## Oticon Opn™ Fitting Range†

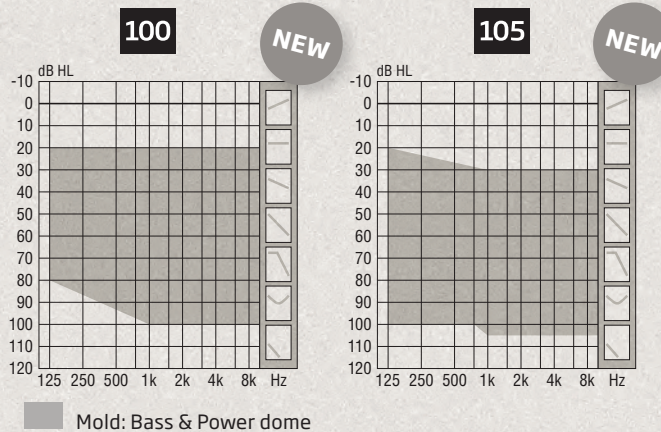

† Fitting range is based on Oticon Opn 1. Details for Oticon Opn 2 & Oticon Opn 3 are available in Technical Data Sheets.

## Small, discreet miniRITE

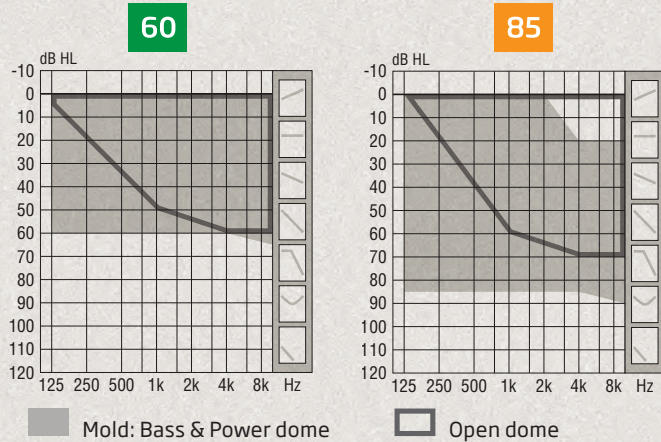

Oticon Opn™ miniRITE has a new, discreet design with a smart single push button for easy operation of volume and programs.

Oticon Opn miniRITE offers clients a discreet hearing aid with a wealth of new features and functionalities including 2.4 GHz wireless technology and Made for iPhone® functionality.\*

Oticon Opn miniRITE uses the proven miniFit receivers and ear pieces, fits up to 105 dB HL and is powered by a 312 battery.

### miniFit receivers and Power Receiver Molds

Select between three receivers with different output performance, with lengths 0-5 for non power speakers and 1-5 for power. The 105 receiver requires a custom Power Receiver Mold.



### Earpieces

miniFit earpieces include standard domes, MicroMold, and LiteTip. Custom earpieces are available in acrylic, soft silicone or VarioTherm® material.



### IP68

Proven IP68 water and dust resistant hearing aids.



	<b>Oticon Opn 1</b>	<b>Oticon Opn 2</b>	<b>Oticon Opn 3</b>
<b>Speech Understanding</b>			
OpenSound Navigator™	Level 1	Level 2	Level 3
Balancing power effect	100%	50%	50%
Max. noise removal	9 dB	5 dB	3 dB
Speech Guard™ LX	Level 1	Level 2	Level 3
Spatial Sound™ LX	4 estimators	2 estimators	2 estimators
Soft Speech Booster LX	•	•	•
<b>Sound Quality</b>			
Clear Dynamics	•	•	-
Binaural Noise Management	•	•	-
Fitting Bandwidth (accessed in software)	10 kHz	8 kHz	8 kHz
Processing Channels	64	48	48
Bass Boost (streaming)	•	•	•
<b>Listening Comfort</b>			
Transient Noise Management	4 configurations	On/Off	On/Off
Feedback shield LX	•	•	•
Wind Noise Management	•	•	•
Binaural Coordination	•	•	•
<b>Personalization &amp; Optimizing Fitting</b>			
YouMatic™ LX	3 configurations	2 configurations	1 configuration
Fitting Bands	16	14	12
Multiple Directionality Options	•	•	•
Adaptation Management	•	•	•
Oticon Firmware Updater	•	•	•
Fitting Formulas	VAC+, NAL-NL1 + 2	VAC+, NAL-NL1 + 2	VAC+, NAL-NL1 + 2
Acoustic Notifications	•	•	•
<b>Connecting to the World</b>			
Stereo streaming (2.4 GHz)	•	•	•
Made for iPhone®	•	•	•
Oticon ON App	•	•	•
ConnectClip (coming in 2017)	•	•	•
Remote Control 3.0	•	•	•
TV Adapter 3.0	•	•	•
AutoPhone	•	•	•

\*Oticon Opn™ is compatible with iPhone 7 Plus, iPhone 7, iPhone SE, iPhone 6s Plus, iPhone 6s, iPhone 6, iPhone 5s, iPhone5c, iPhone 5, 9.7-inch iPad Pro, 12.9-inch iPad Pro, iPad Air 2, iPad Air (4<sup>th</sup> generation), iPad mini 4, iPad mini 3, iPad mini 2, iPad mini, and iPod touch 5<sup>th</sup> and 6<sup>th</sup> generation. Devices must be running iOS 9.3 or later. Please visit [www.oticon.com](http://www.oticon.com) for more details.

## OTICON | Opn

# Product comparison



 **TwinLink**  
NFMI + 2.4 GHz

**oticon**  
PEOPLE FIRST